



中国认可  
国际互认  
检测  
TESTING  
CNAS L1842



佛山中纺联检验技术服务有限公司  
CNTAC Testing Service Co.,Ltd.(Foshan)

<b>Testing Report</b>	Security website: <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a> Security code: 0150221143
Report No: ZFLJ2640106	Page 1 of 3



### Applicant Information

Applicant Name : SANKOM SA  
Applicant Address : SANKOM Building Rue de la Moliere 2 CH-2800 Delemont

### Sample Information

Sample Description : AVEVITTA PROTECT 2.0  
Brand : AVEVITTA  
Sample Quantity : 20 pieces  
Model : AV-020

- Sample Receiving Date : 2020-08-01  
- Report Date : 2020-08-06  
- The original sample is stucked on the last paper.

### Test Performed

Judgement according to:  
EN 14683:2019+AC:2019 Medical face masks – Requirements and test methods Type II R  
- Selected test(s) as requested by applicants. For details, refer to attached page(s).


### Pronounce

The results shown in this report refer only to sample(s) tested unless otherwise stated.  
All the items are performed in the standard conditions, except the noted cases.  
Except for the requirement of the client, the test results and the conformity judgement of this report do not take the uncertainty of the test results into account.

Signed for and on behalf of  
CNTAC Testing Service Co.,Ltd.(Foshan)

Approved by

罗拉莲

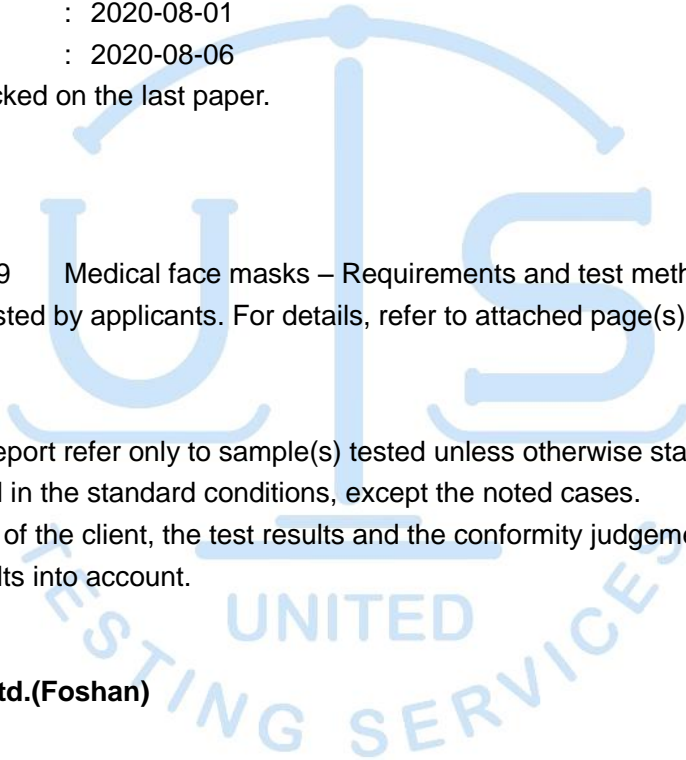



中国认可  
国际互认  
检测  
TESTING  
CNAS L1842



<b>Testing Report</b>	Security website: <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a> Security code: 0150221143
Report No: ZFLJ2640106	Page 2 of 3



--- Test Result ---

Test Result	Requirements	Judgement	
<b>1. Bacterial Filtration Efficiency</b> <u>EN 14683:2019+AC: 2019 Appendix B</u>			
Flow rate: 28.3L/min, Test area: 77cm <sup>2</sup> , Test bacteria culture: Staphylococcus aureus ATCC 6538, Inside of the test mask was facing towards the challenge aerosol, The average plate count results of the positive controls: 2.5×10 <sup>3</sup> CFU, The average plate count results of the negative controls: <1CFU			
Unit: <%>			
sample 1#	99.8	≥ 98	Pass
sample 2#	99.9		
sample 3#	99.8		
sample 4#	99.8		
sample 5#	99.9		
Minimum value	99.8		





中国认可  
国际互认  
检测  
TESTING  
CNAS L1842



佛山中纺联检验技术服务有限公司  
CNTAC Testing Service Co.,Ltd.(Foshan)

<b>Testing Report</b>	Security website: <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a> Security code: 0150221143
Report No: ZFLJ2640106	Page 3 of 3



**Original Sample**



===== End of Report =====

