



Test Report

No.: T32020260713SN

Date: JUL 20, 2020

Page 1 of 3

SANKOM SA
SANKOM BUILDING RUE DE LA MOLIERE 2, DELEMONT, CH-2800, SWITZERLAND

The following samples were submitted and identified on behalf of the client as:

AVEVITTA PROTECT 2.0

SGS Case No. : CA320202639699
 Lot No. / Batch Code : NOT PROVIDED
 Sample Description : FABRIC MASK IN ASSORT COLOUR
 Quantity Submitted : 22 PCS
 Style / Item No. : AV-020
 Manufacturer : SANKOM SA
 Country of Origin : CHINA

Sample Receiving Date : JUN 29, 2020
 Test Performing Date : JUN 29, 2020 TO JUL 20, 2020

Test Requested : Please refer to the result summary.

Test Method & Results : Please refer to next page(s).

Result Summary :

Test Requested	Result Summary
Determination of Breathability (EN14683:2019+AC:2019 Appendix C Differential pressure) (For Blue Style only)	See Result

Signed for and on behalf of
SGS Hong Kong Ltd.

Au Kam Chi, Gigi
Technical Manager

This document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Documents.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law.

Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.





Test Report

No.: T32020260713SN

Date: JUL 20, 2020

Page 2 of 3

Determination of Breathability (EN14683:2019+AC:2019 Appendix C Differential pressure)


Test Side : Blue Colour Without Print (Inside)
Test Location : Geometry Centre
Pre-Conditioning : Minimum of 4 hours at 21±5°C and 85±5% R.H.
Test Area : 4.9 cm²
Flow Rate : 8 l/min

Test Specimen	ΔP (mm H ₂ O/cm ²)	ΔP (Pa/cm ²)
1	2.9	28.3
2	2.4	23.9
3	3.5	34.6
4	3.0	29.2
5	3.5	34.4
Average	3.06	30.08

This document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Documents.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law.

Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.

Photo Appendix

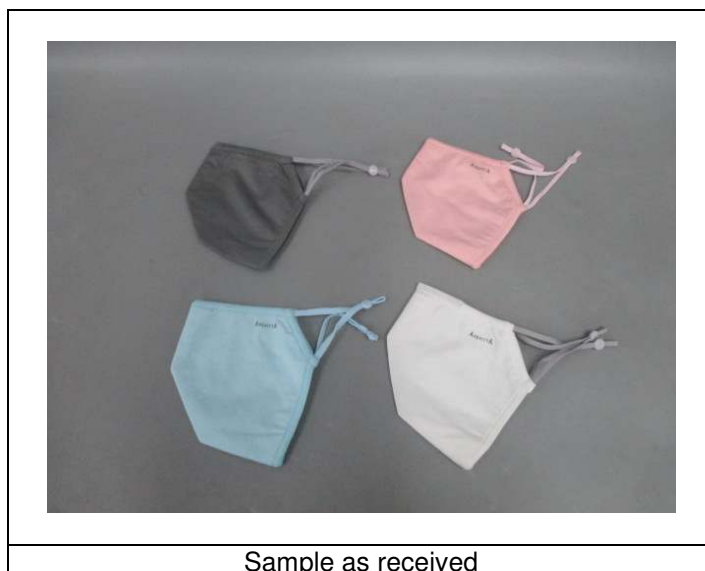

*** End of Report ***

This document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at <http://www.sgs.com/en/Terms-and-Conditions.aspx> and, for electronic format documents, subject to Terms and Conditions for Electronic Documents at <http://www.sgs.com/en/Terms-and-Conditions/Terms-e-Documents.aspx>. Attention is drawn to the limitation of liability, indemnification and jurisdiction issues defined therein. Any holder of this document is advised that information contained hereon reflects the Company's findings at the time of its intervention only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not exonerate parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law.

Unless otherwise stated the results shown in this test report refer only to the sample(s) tested and such sample(s) are retained for 30 days only.