



中国认可  
国际互认  
检测  
TESTING  
CNAS L1842



佛山中纺联检验技术服务有限公司  
CNTAC Testing Service Co.,Ltd.(Foshan)

<b>Testing Report</b>	Security website: <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a> Security code: 4303853843
Report No: ZFLJ2644809	Page 1 of 3



### Applicant Information

Applicant Name : SANKOM SA  
Applicant Address : SANKOM Building Rue de la Moliere 2 CH-2800 Delemont

### Sample Information

Sample Description : AVEVITTA PROTECT 2.0  
Brand : AVEVITTA  
Sample Quantity : 34 pieces  
Model : AV-020

- Sample Receiving Date : 2020-08-28  
- Report Date : 2020-08-29  
- The original sample is stucked on the last paper.

### Test Performed

Judgement according to:  
ASTM F2100-2019e1 Standard Specification for Performance of Materials Used in Medical Face Masks  
Level 1  
- Selected test(s) as requested by applicants. For details, refer to attached page(s).


### Pronounce

The results shown in this report refer only to sample(s) tested unless otherwise stated.  
All the items are performed in the standard conditions, except the noted cases.  
Except for the requirement of the client, the test results and the conformity judgement of this report do not take the uncertainty of the test results into account.

Signed for and on behalf of  
CNTAC Testing Service Co.,Ltd.(Foshan)

Approved by

张志荣

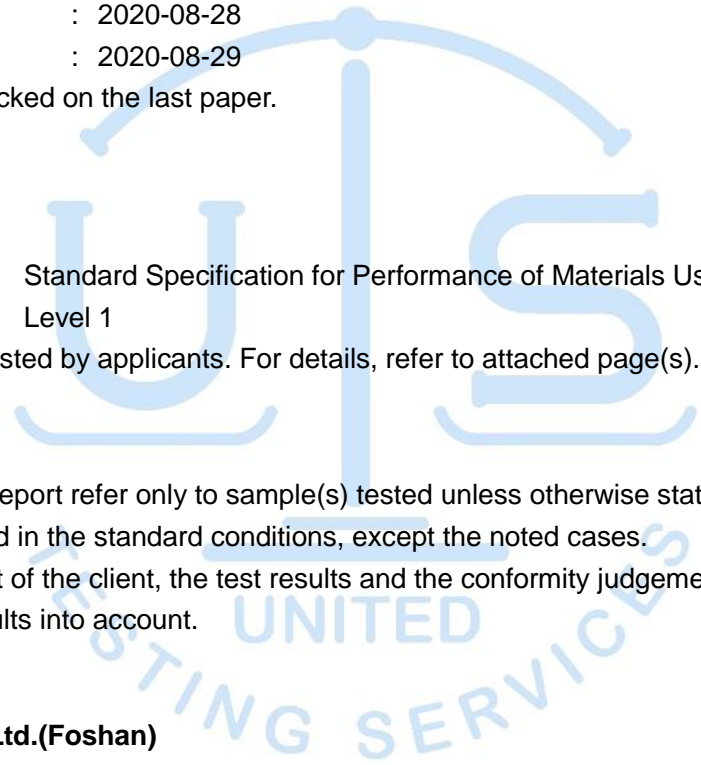



中国认可  
国际互认  
检测  
TESTING  
CNAS L1842



佛山中纺联检验技术服务有限公司  
CNTAC Testing Service Co.,Ltd.(Foshan)

<b>Testing Report</b>	Security website: <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a> Security code: 4303853843
<b>Report No: ZFLJ2644809</b>	Page 2 of 3



--- Test Result ---

Test Result	Requirements	Judgement
<b>1. Resistance to Penetration by Synthetic Blood</b> <i>ASTM F1862/F1862M-17</i>		
Synthetic blood surface tension: 0.040N/m, Distance between blow head front end and target area: 300mm, Injection synthetic blood volumes: 2mL, Test Blood Pressure: 80mmHg, Injection velocity: 450cm/s, Use a fixed target board .		
sample 1#	None seen	Pass Pressure at 80mmHg
sample 2#	None seen	Pass
sample 3#	None seen	
sample 4#	None seen	
sample 5#	None seen	
sample 6#	None seen	
sample 7#	None seen	
sample 8#	None seen	
sample 9#	None seen	
sample 10#	None seen	
sample 11#	None seen	
sample 12#	None seen	
sample 13#	None seen	
sample 14#	None seen	
sample 15#	None seen	
sample 16#	None seen	
sample 17#	None seen	
sample 18#	None seen	
sample 19#	None seen	
sample 20#	None seen	
sample 21#	None seen	
sample 22#	None seen	
sample 23#	None seen	
sample 24#	None seen	
sample 25#	None seen	
sample 26#	None seen	
sample 27#	None seen	
sample 28#	None seen	
sample 29#	None seen	
sample 30#	None seen	
sample 31#	None seen	
sample 32#	None seen	





中国认可  
国际互认  
检测  
TESTING  
CNAS L1842



佛山中纺联检验技术服务有限公司  
CNTAC Testing Service Co.,Ltd.(Foshan)

<b>Testing Report</b>	Security website: <a href="http://www.fcl-sz.org.cn">www.fcl-sz.org.cn</a> Security code: 4303853843
Report No: ZFLJ2644809	Page 3 of 3



**Original Sample**



=====**End of Report**=====

